

THE SILENTIA SCREEN SYSTEM



# The leading alternative for optimal hygiene and patient privacy

CE





We set the standard in 1990.  
Today we have delivered 100 000  
screens around the world.

Silentia offers a complete solution to meet screening requirements in healthcare. Because the individual screens in the range are compatible with each other and have the EasyClick feature, you can easily move a folding screen from a wall to a trolley – or combine a fixed and folding screen to get maximum benefit from all the system possibilities.

#### **EasyClean™ for optimal hygiene**

Screens have smooth surfaces that can withstand the disinfectants used in care situations today. The enclosed wheels can easily be removed and washed in a dishwasher. All products meet strict hygiene requirements and all surfaces are accessible for cleaning – that's EasyClean.

#### **Folds back easily with one hand – EasyReturn™**

Patented hinges have a *Built-in Memory* which enables easy handling and prevention of cross infection due to having only one contact spot: its ergonomically designed handle.

#### **Fixed or mobile – the choice is yours with EasyClick™**

Moveability is one click away. A wall-mounted folding screen can be mounted on a trolley for mobile use, or on other system components. There are many possibilities, making this a complete system, easy to use in all departments.


#### **CE mark in accordance to MDD 93/42/EEC**

Silentia products meets the essential requirements in Annex 1 of the Medical Devices Directive (MDD).



- Several lengths and heights are available
- Extreme lightweight construction
- Made from environmentally friendly materials
- Section panels are made of high pressure laminate or polycarbonate profiles which withstand the disinfectants that are used within the hospital.


**EasyClean**

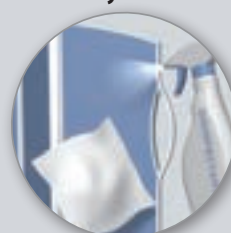


**EasyClean**



**EasyReturn**

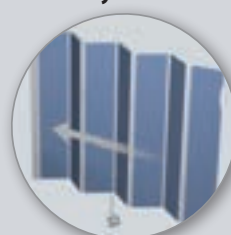

**EasyClick**

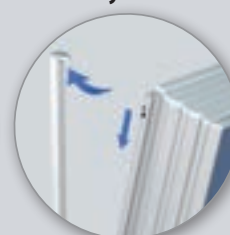




## Folding screens

Simplicity is one of our guiding principles in product development. It has to be easy to screen off patients when they need privacy or simply want time to themselves. It also has to be easy to change things around and expand the system – so all components are designed to be compatible with each other. See pages 4–5.

### STANDARD VERSION – MOUNTED ON A WALL

A wall mounting permits easy fitting. EasyClick screen mounting.

The screen occupies only 26 cm of wall space.

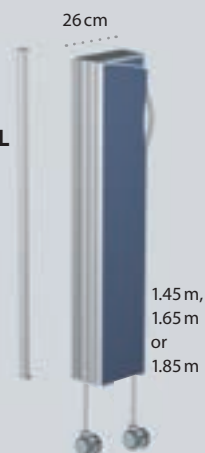

### SCREEN TROLLEY

Accessory for folding screens. EasyClick screen mounting.

The wheel base takes a minimum of space.




### POLE FOR FLOOR MOUNTING

Designed to be anchored to the floor. This is applied when a free-standing screen is required. EasyClick screen mounting.




### WALL-MOUNTED WITH SPACER

Used where a screen has to be mounted a short distance from the wall, to provide clearance for a care panel, for example. EasyClick screen mounting.

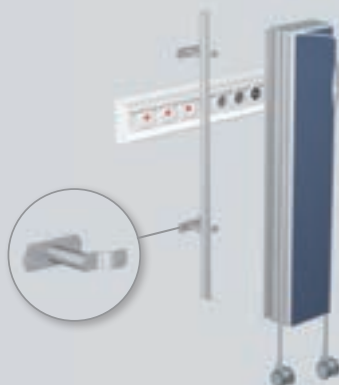

### MOUNTED ON WALL RAIL

A clamp is used to mount a screen on a wall rail so that it can be adjusted sideways. Two wall rails should be used if the screen is longer than 2.25 m. EasyClick screen mounting.

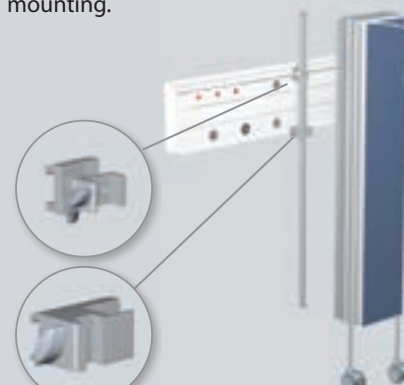







# Mobile folding screens

## MOUNTED ON TROLLEY

The screen can be mounted on a trolley to provide convenient mobile screening. This trolley is suitable for folding screens of various lengths: from 5 to 15 sections. Height: 1.45 m, 1.65 m or 1.85 m. The patented design takes a minimum of floor space while providing maximum stability.

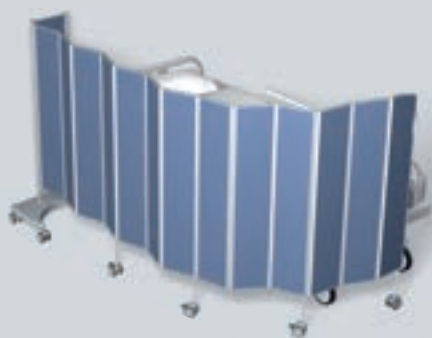
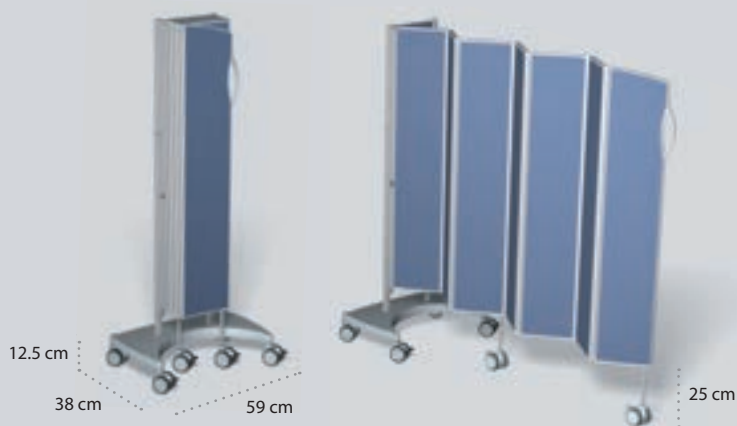




## Mobile bed end screens



### SILENTIA BED END SCREEN 1.5 M

Provides effective screening at the end of the bed – or for changing. Extended width: 1.5 m. Height: 1.55 m or 1.85 m. The patented design takes a minimum of floor space while providing maximum stability.

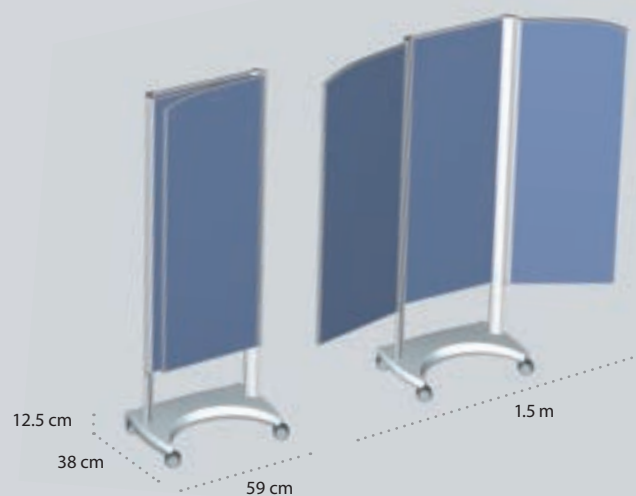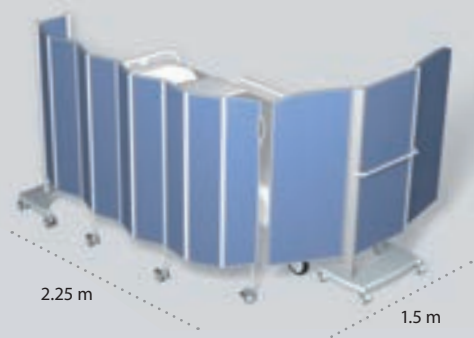

### SILENTIA BED END SCREEN 2.0 M

This new bed end screen is ideal when longer screening is required. Extended width: 2.0 m. Height: 1.55 m or 1.85 m.

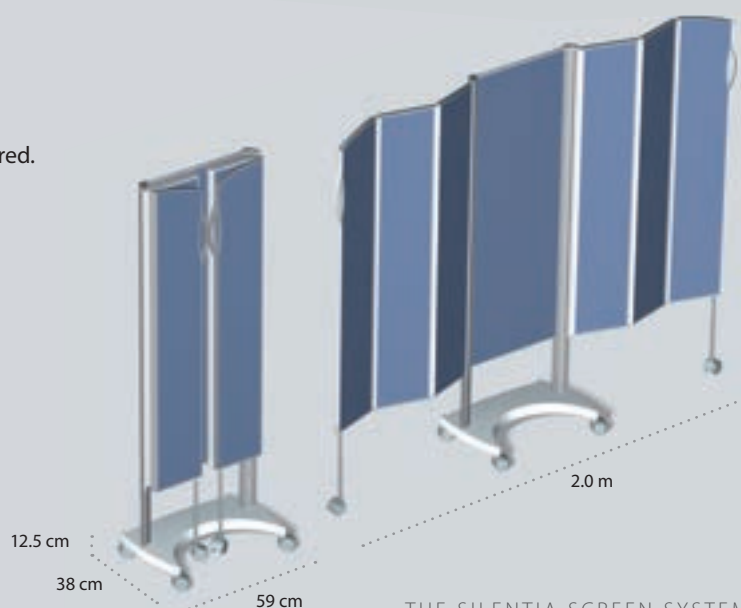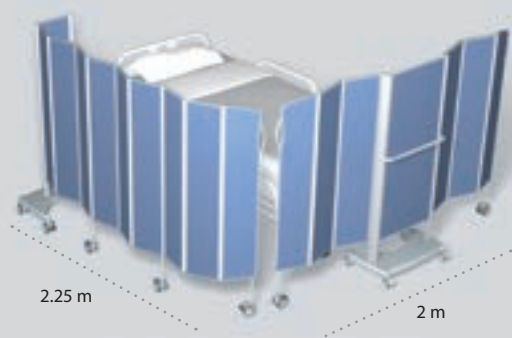



## Fixed screens


Our fixed screens have been developed for applications where more permanent screening is required. This system also permits a wider range of possibilities. A fixed screen between beds can, for example, be combined with a pull-out screen at the foot of the bed. Fixed screens can be used to create changing cubicles and other types of room within a room, without making alterations to walls.

### STANDARD SIZES

LENGTH	HEIGHT	WEIGHT	ART. NO. COLOURED	ART. NO. DAYLIGHT TOP	ART. NO. DAYLIGHT TOP DUO	ART. NO. DAYLIGHT	ART. NO. DAYLIGHT DUO
1.3 m	1.45 m	6.2 kg	6150	6150-DT	—	6150-D	—
1.3 m	1.65 m	7.2 kg	6250	6250-DT	—	6250-D	—
1.3 m	1.85 m	8.3 kg	6350	6350-DT	—	6350-D	—
2.4 m	1.45 m	11.4 kg	6155	6155-DT	6155-DTD	6155-D	6155-DD
2.4 m	1.65 m	13.3 kg	6255	6255-DT	6255-DTD	6255-D	6255-DD
2.4 m	1.85 m	15.3 kg	6355	6355-DT	6355-DTD	6355-D	6355-DD



DAYLIGHT, FIXED

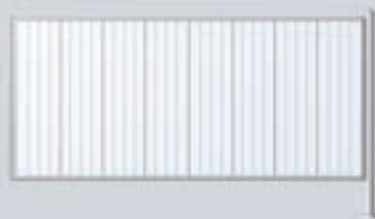


DAYLIGHT TOP DUO, FIXED



DAYLIGHT DUO, FIXED



COLOURED SCREEN, FIXED

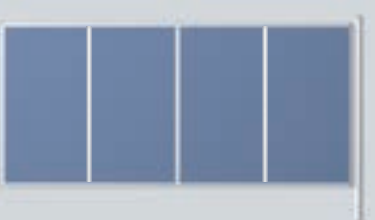

DAYLIGHT TOP, FIXED







With floor poles you can create a service corridor behind the bed head



## SILENTIA FLOOR MOUNTING

A standard floor pole (**A**) is used to install a single fixed screen with EasyClick. If you want the pole to serve as a mounting point for additional fixed or folding screens there are special versions (**B-D**) with up to four mounting points.

The pole is normally secured to the floor for maximum rigidity – this is especially important if the pole supports more than one screen. If the pole is only intended to support a single wall-mounted screen it can simply rest on the floor.

## STANDARD SIZES FLOOR POLES

HEIGHT	WEIGHT	ART. NO.	ART. NO.	ART. NO.	ART. NO.
1.45 m	3.3 kg	6561-A	6561-B	6561-C	6561-D
1.65 m	3.5 kg	6562-A	6562-B	6562-C	6562-D
1.85 m	3.8 kg	6563-A	6563-B	6563-C	6563-D

Floor pole attachment possibilities – seen from above.



The pole and the cover of the pole at the base are attractively designed with a smooth surface that is easy to clean.




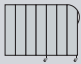




## EXTENSION SCREEN

Moves the mounting point for the folding screen away from the wall, either to make space for the bedside table, or to extend a folding screen. Available in four lengths and three heights.

HEIGHT	ART. NO. LENGTH 0.5 M	ART. NO. LENGTH 0.75 M	ART. NO. LENGTH 1.0 M	ART. NO. LENGTH 1.2 M
1.45 m	6162	6163	6164	6165
1.65 m	6262	6263	6264	6265
1.85 m	6362	6363	6364	6365

# Silentia folding screens

## STANDARD SIZES

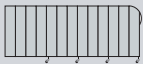	NO. OF PANELS	ART. NO. COLOURED	ART. NO. DAYLIGHT TOP	ART. NO. DAYLIGHT DUO	ART. NO. DAYLIGHT	LENGTH	HEIGHT	WEIGHT	FOLDED SIZE	
									WIDTH	DEPTH
	5	6105	6105-DT	6105-DD	6105-D	1.25 m	1.45 m	8 kg	0.25 m	0.1 m
		6205	6205-DT	6205-DD	6205-D	1.25 m	1.65 m	9.5 kg	0.25 m	0.1 m
		6305	6305-DT	6305-DD	6305-D	1.25 m	1.85 m	11.5 kg	0.25 m	0.1 m
	7	6107	6107-DT	6107-DD	6107-D	1.75 m	1.45 m	12 kg	0.25 m	0.14 m
		6207	6207-DT	6207-DD	6207-D	1.75 m	1.65 m	14 kg	0.25 m	0.14 m
		6307	6307-DT	6307-DD	6307-D	1.75 m	1.85 m	16 kg	0.25 m	0.14 m
	9	6109	6109-DT	6109-DD	6109-D	2.25 m	1.45 m	15.5 kg	0.25 m	0.18 m
		6209	6209-DT	6209-DD	6209-D	2.25 m	1.65 m	18 kg	0.25 m	0.18 m
		6309	6309-DT	6309-DD	6309-D	2.25 m	1.85 m	20.5 kg	0.25 m	0.18 m
	11	6111	6111-DT	6111-DD	6111-D	2.75 m	1.45 m	18.5 kg	0.25 m	0.22 m
		6211	6211-DT	6211-DD	6211-D	2.75 m	1.65 m	21.5 kg	0.25 m	0.22 m
		6311	6311-DT	6311-DD	6311-D	2.75 m	1.85 m	24.5 kg	0.25 m	0.22 m
	13	6113	6113-DT	6113-DD	Not available	3.25 m	1.45 m	22 kg	0.25 m	0.26 m
		6213	6213-DT	6213-DD		3.25 m	1.65 m	25.5 kg	0.25 m	0.26 m
		6313	6313-DT	6313-DD		3.25 m	1.85 m	29 kg	0.25 m	0.26 m
	15	6115	6115-DT	6115-DD	Not available	3.75 m	1.45 m	25.5 kg	0.25 m	0.3 m
		6215	6215-DT	6215-DD		3.75 m	1.65 m	29.5 kg	0.25 m	0.3 m
		6315	6315-DT	6315-DD		3.75 m	1.85 m	33.5 kg	0.25 m	0.3 m


The table shows our standard sizes in 0.5 m length increments. Screens are also available in increments of 0.25 m if the total length is between 0.75 m and 3.75 m. The standard heights – 1.45 m, 1.65 m and 1.85 m – can be lowered by – 0.1 m by determining the height of the wheel legs. The height is decided when you order your Silentia folding screen.

**Note:** The longest screens (12–15 sections) are intended for use in an L-shape around the bed.



## CHOOSE THE RIGHT MOUNTING TO SUIT YOUR SCREENING NEED

(A) The folding screen mounting is normally on the left side, which means that the screen is folded to the right of the wall mounting. (B) To prevent the screen from taking up too much space next to a doorway, for example, screens can also be ordered with the mounting on the opposite side.

A

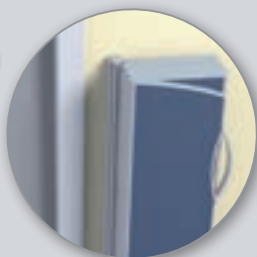

B

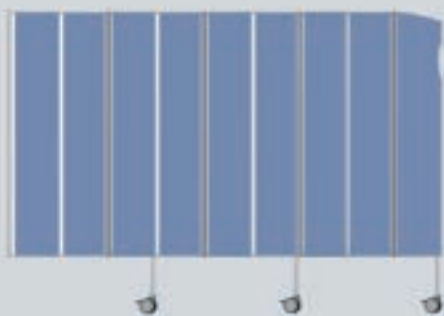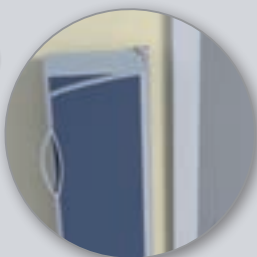

**COLOURED SCREENS**



**DAYLIGHT TOP**

Top section is about 0.3–0.4 m high.



**DAYLIGHT DUO**


Frosted sections are in pairs.



**DAYLIGHT**



# Colors and patterns

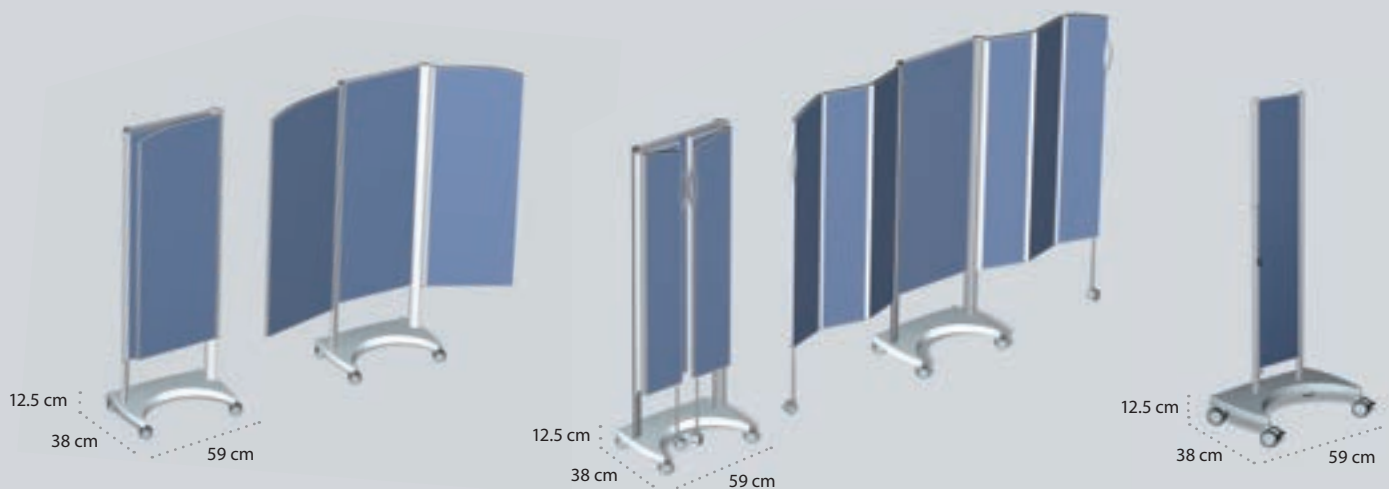

## SILENTIA BED END SCREEN 1.5 M

Art. no. 6522, height 1.55 m

Art. no. 6523, height 1.85 m

Extended width 1.5 m, folded width 0.55 m.

This screen is intended for situations where a free-standing screen is needed, e.g. for changing or at the end of a bed. The wheeled stand is low in height and easy to use.

## SILENTIA BED END SCREEN 2.0 M

Art. no. 6525, height 1.55 m

Art. no. 6526, height 1.85 m

Extended width 2.0 m, folded width 0.55 m.

Provides effective privacy screening at end of bed, in corridors or for changing. The wheels can be locked when the screen is extended.

The wheeled stand is low in height and designed for optimum performance.

## SCREEN TROLLEY

Art. no. 6514, height 1.45 m

Art. no. 6518, height 1.65 m

Art. no. 6520, height 1.85 m

Accessory for folding screens up to 15 sections long.

EasyClick screen mounting.

**Note:** The longest screens (12–15 sections) are intended for use in an L-shape around the bed.

# Clamps and spacer



## STANDARD CLAMP FOR WALL RAIL

Art. no. 6502

For attaching to wall rail. A flexible solution that enables the screen to be moved sideways along the wall rail. Distance between wall rail and screen: 40 mm.



## SPACER 20 MM

Art. no. 6506

Used to move a screen away from the wall, to make space for a moulding, for example. Supplied in pairs.



## SPACER 100 OR 200 MM

Art. no. 6507, 100 mm

Art. no. 6508, 200 mm

Used to move a screen away from the wall. Supplied in pairs.



## ADJUSTABLE SPACER

Art. no. 6509, 200–350 mm

Art. no. 6510, 300–500 mm

Used to move a screen well away from the wall. Supplied individually.



## SUPPORT CLAMP

Art. no. 6524

Used with standard clamp 6502 where two wall rails are mounted on the wall.



## CLAMP WITH SPACER

Art. no. 6535

A useful solution where space is required between a screen and a wall rail, such as in dialysis departments where panel access is required. Distance between wall rail and screen: 150 mm.



## ADJUSTABLE CLAMP

Art. no. 6516

Adjusts to suit variable spacing between upper and lower wall rails. Clamps can be adjusted between 0 and 150 mm. Used together with standard clamp or clamp with spacer.



## BACK SUPPORT ON WALL

Art. no. 6519

Used as backward support when using a clamp mounted on a wall rail.

# Wall mounting and floor pole



## WALL MOUNTING

Art. no. 6539, height 1.45 m

Art. no. 6544, height 1.65 m

Art. no. 6540, height 1.85 m

For convenient mounting direct on the wall. One wall mounting is normally supplied with the screen. Additional wall mountings can be ordered separately.



## WALL MOUNTING WITH FLOOR SUPPORT

Art. no. 6537, height 1.45 m

Art. no. 6545, height 1.65 m

Art. no. 6538, height 1.85 m

Wall mounting that provides support from the floor when fitted with a clamp.



## WALL MOUNTING WITH FLOOR SUPPORT AND TOP EXTENSION

Art. no. 6541, height 1.45 m

Art. no. 6546, height 1.65 m

Art. no. 6542, height 1.85 m

Wall mounting used with clamp where a wall rail crosses above the screen.



## FLOOR POLE

Art. no. 6561A, height 1.45 m

Art. no. 6562A, height 1.65 m

Art. no. 6563A, height 1.85 m

The floor pole can be used to support a wall-mounted screen or as a mounting point for folding screens. The pole is screwed securely to the floor. Other versions can be seen on page 7.

# Folding screens on wall rail




## SUPPORT CLAMP (ABOVE)

Art. no. 6524


## STANDARD CLAMP (BELOW)

Art. no. 6502

This support clamp combination is chosen when 2 rails are in a vertical line.

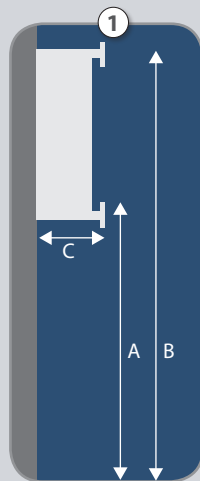

IMPORTANT AND NECESSARY MEASUREMENTS

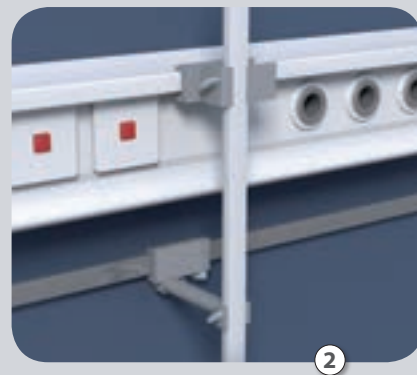


## STANDARD CLAMP (ABOVE)

Art. no. 6502

## ADJUSTABLE CLAMP (BELOW)

Art. no. 6516

This support clamp combination is chosen when you have 1 rail closer to the wall than the other.

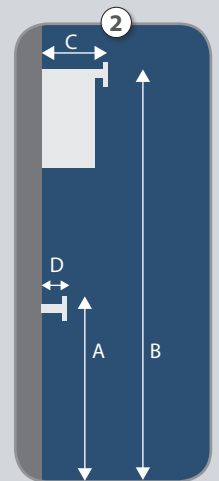

IMPORTANT AND NECESSARY MEASUREMENTS



## CLAMP WITH SPACER (ABOVE AND BELOW)

Art. no. 6535

This combination is needed when you need a space between the rail and the screen for tubes and other equipment. Important measurements see no. 2.



1,85 cm




1,45 cm



## STANDARD CLAMP AND BACK

SUPPORT ON WALL Art. no. 6502 + 6519

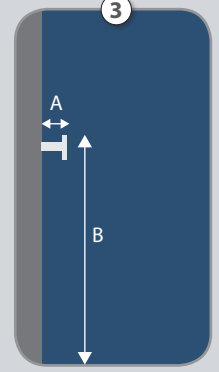

IMPORTANT AND NECESSARY MEASUREMENTS

# Wall mount spacer



## SPACER

100 mm – Art. no. 6507

200 mm – Art. no. 6508

Moves the screen 100 mm or 200 mm from the wall. Suitable for use when a supply panel, electric socket, gas point or switch is located behind the screen. Also provides space for hoses and cables behind the screen.




## ADJUSTABLE SPACER


200 – 350 mm – Art. no. 6509

300 – 500 mm – Art. no. 6510

Flexible spacer for use when a screen needs to be moved well away from the wall, for example if a shelf or deep supply panel is located behind the screen.

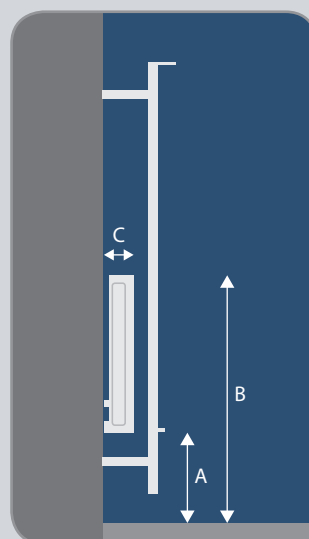

IMPORTANT AND NECESSARY MEASUREMENTS FOR RADIATORS (A, B, AND C)

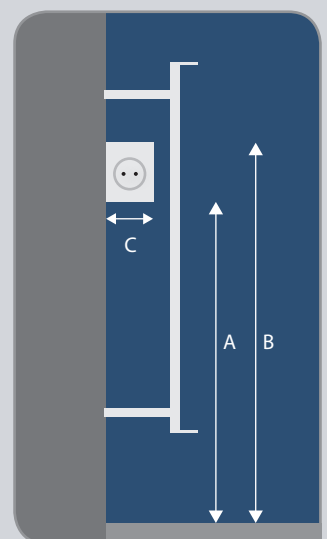

IMPORTANT AND NECESSARY MEASUREMENTS FOR ELECTRICAL OUTLETS (A, B, AND C)





CE

**Ansvarig utgivare** FOTO: ANDERS SÄLLSTRÖM. TRYCK: STAR-TRYCK 11.14